

## I/O Ports

The Medical-Grade Integrated LCD System has a comprehensive set of I/O ports. The following ports are located along the bottom side of the unit.


Front




Back



Bottom

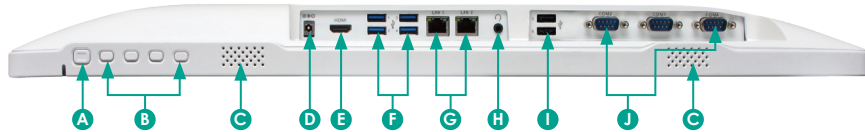

|                               |  |
|-------------------------------|--|
| <b>A</b> Power button         | <b>F</b> USB 3.0 ports                               |
| <b>B</b> Programmable buttons | <b>G</b> Ethernet ports (RJ45)                       |
| <b>C</b> Speakers             | <b>H</b> Audio jack                                  |
| <b>D</b> DC-in                | <b>I</b> USB 2.0 ports                               |
| <b>E</b> HDMI output          | <b>J</b> COM ports                                   |
|                               | <b>K</b> Potential equalization conductor (optional) |

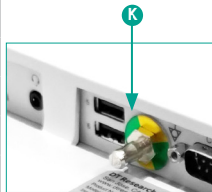

### Cleaning the Screen

- A screen cleaning solution that is alcohol-free and non-abrasive can be used to clean the touch screen.
- Cleaning with a micro-fiber cloth is recommended.
- Please spread the solution onto the cloth and then clean the touch screen.

### Cleaning the Anti-microbial Enclosure

- Use a soft/non-abrasive cloth moistened with water to clean the enclosure.
- If using a cleaner, an alcohol-free and oxide-free cleaning liquid is recommended.
- To prevent scratching the anti-microbial coating, please wipe gently.



MEDICAL – GENERAL MEDICAL EQUIPMENT  
AS TO ELECTRICAL SHOCK, FIRE AND MECHANICAL  
HAZARDS ONLY IN ACCORDANCE WITH ANSI/AAMI  
ES60601-1 (2005) + AMD (2012) & "CAN/CSA-C22.2  
No. 6060-1 (2008) + (2014)



DT Research, Inc.  
2000 Concourse Drive, San Jose, CA 95131  
Copyright © 2019, DT Research, Inc. All Rights Reserved.  
BOG031319 DT524T ENG

# DT Research

# Medical-Grade Integrated LCD System

## DT524T

### BASIC OPERATION GUIDE



### ENGLISH

### INTRODUCTION

Thank you for acquiring DT Research's Medical-Grade Integrated LCD System. With a 24-inch TFT display and powered by the Intel® Core™ i processor, the Medical-Grade Integrated LCD System offers optimal combinations of performance and power savings. With fully-integrated point-of-care modules within an elegant, space-saving design, the Medical-Grade Integrated LCD System is the optimum solution to enhance workflow and service.

Please take a few moments to review the contents of this document to ensure that the setup and startup proceed smoothly. The Medical-Grade Integrated LCD System is ready for use, out of the box, in its default configuration when powered by the power source provided. The following documentation offers guidance on the hardware elements and features of the computer. Please refer to your device provider for information pertaining to the software operating system or software applications.

### Package Contents

- DT524T Medical-Grade Integrated LCD System
- AC-DC power adapter with power cord
- Basic operation guide

### Major Options

The Medical-Grade Integrated LCD System has a number of optional modules. These include:

- 2 megapixel front camera
- UPS battery
- Desktop stand

## Precautions

- Always exercise care when operating and handling the Medical-Grade Integrated LCD System.
- Never disassemble any portion of the enclosure, as this will void any product warranty on the Medical-Grade Integrated LCD System.
- Do not use any AC/DC adapter other than the one provided with the device or a replacement acquired from the manufacturer.
- In the unlikely event that smoke, abnormal noise or strange odor is present, immediately power down the Medical-Grade Integrated LCD System and disconnect all power sources. Please report the problem to your device provider immediately.

## Basic Features

The Medical-Grade Integrated LCD System integrates a bright 24" display with a high performance system, USB ports, and integrated options such as camera, RFID reader and UPS battery for a comprehensive point-of-care solution.

### Attach The Tabletop Stand (optional)

1. Insert the top two tabs of the tabletop stand into the top two open slots on the DT524T back panel.
2. Push the latch to pivot the two lower tabs of the tabletop stand.
3. Push the two lower tabs into the two lower slots on the DT524T back panel. Make sure the latch locks back into position so the stand is securely attached.



### Detach the Tabletop Stand (optional)

1. Push the latch to pivot the two lower tabs.
2. Pull out the two lower tabs, then the two upper tabs from the DT524T back panel.



## Powering ON and OFF

Please use the AC-DC adapter with the Medical-Grade Integrated LCD System for the power supply. To activate the Medical-Grade Integrated LCD System, push and quickly release the Power Button and the display will come on in a few seconds. To put in Standby mode, push and quickly release the Power Button. To turn off for extended storage, power off the device safely using any software function that “shuts down computer” provided in the software operating system.

### NOTE:

To conserve power, use (push and quick release) the Power Button to put the device in “Standby” mode while not in use. Pushing briefly on the same button will wake up the system within seconds.

### NOTE:

Avoid using the Power Button (“hold 4+ seconds” feature) to turn off the device—this form of hardware shutdown is intended to be a means of recovery from device lockups, and not as normal operation.

## Wireless Networking

### Wireless LAN

The Medical-Grade Integrated LCD System is often delivered with an embedded (user-inaccessible) 802.11ac WLAN adapter equipped with a hidden custom antenna.

- Through the support of typical WLAN adapters, the Medical-Grade Integrated LCD System should be able to detect all 802.11 access points in the vicinity for you to select the access point of your choice for connection.
- The SSID and WEP/WPA/WPA2 (if enabled) parameters on the Medical-Grade Integrated LCD System and the access points have to match. The SSID is case-sensitive and it is recommended that you enable WEP/WPA/WPA2 encryption (or advanced alternatives) for secure access.
- When WEP/WPA/WPA2 is enabled, you may need to consult your network administrator or your networking equipment literature to properly configure associated settings such as Authentication mode, etc.
- Refer to the access point operating manuals for setting up the 802.11 access points.

### Bluetooth

The Medical-Grade Integrated LCD System features a built-in Bluetooth adapter that operates on the Microsoft Windows Bluetooth protocol. The Bluetooth configuration application is invoked from the System Tray or from the Control Panel. Follow the instructions and options provided within the application to configure and invoke Bluetooth connectivity with the corresponding peripherals.

### NOTE:

Bluetooth devices or accessories that are not compatible with the Microsoft Windows Bluetooth protocol may not work with the Medical-Grade Integrated LCD System.